



This is a digital copy of a book that was preserved for generations on library shelves before it was carefully scanned by Google as part of a project to make the world's books discoverable online.

It has survived long enough for the copyright to expire and the book to enter the public domain. A public domain book is one that was never subject to copyright or whose legal copyright term has expired. Whether a book is in the public domain may vary country to country. Public domain books are our gateways to the past, representing a wealth of history, culture and knowledge that's often difficult to discover.

Marks, notations and other marginalia present in the original volume will appear in this file - a reminder of this book's long journey from the publisher to a library and finally to you.

Usage guidelines

Google is proud to partner with libraries to digitize public domain materials and make them widely accessible. Public domain books belong to the public and we are merely their custodians. Nevertheless, this work is expensive, so in order to keep providing this resource, we have taken steps to prevent abuse by commercial parties, including placing technical restrictions on automated querying.

We also ask that you:

- + *Make non-commercial use of the files* We designed Google Book Search for use by individuals, and we request that you use these files for personal, non-commercial purposes.
- + *Refrain from automated querying* Do not send automated queries of any sort to Google's system: If you are conducting research on machine translation, optical character recognition or other areas where access to a large amount of text is helpful, please contact us. We encourage the use of public domain materials for these purposes and may be able to help.
- + *Maintain attribution* The Google "watermark" you see on each file is essential for informing people about this project and helping them find additional materials through Google Book Search. Please do not remove it.
- + *Keep it legal* Whatever your use, remember that you are responsible for ensuring that what you are doing is legal. Do not assume that just because we believe a book is in the public domain for users in the United States, that the work is also in the public domain for users in other countries. Whether a book is still in copyright varies from country to country, and we can't offer guidance on whether any specific use of any specific book is allowed. Please do not assume that a book's appearance in Google Book Search means it can be used in any manner anywhere in the world. Copyright infringement liability can be quite severe.

About Google Book Search

Google's mission is to organize the world's information and to make it universally accessible and useful. Google Book Search helps readers discover the world's books while helping authors and publishers reach new audiences. You can search through the full text of this book on the web at <http://books.google.com/>

24503385315



LANE MEDICAL LIBRARY STANFORD STOR
0211 .B35 1885
Manipulation of the microscope / by Edwin

MANIPULATION

MICROSCOPE 

EDWARD BAISCH.

0211
B35
1885

✱LIBRARY✱

OF

Cooper Medical College

DATE

NO. 96

SHELF B6

GIFT OF

LANE

MEDICAL



LIBRARY

LEVI COOPER LANE FUND

1910

✧ LIBRARY ✧
OF
Cooper Medical College

DATE

No. 96 SHELF B6

GIFT OF



1930





MANIPULATION
OF THE
MICROSCOPE,

LANE LIBRARY

EDWARD BAUSCH.

ILLUSTRATED.



PUBLISHED BY BAUSCH & LOMB OPTICAL CO.

ROCHESTER, N. Y. :
POST-EXPRESS BOOK AND JOB PRINTING HOUSE.
1885.

At



YASUJI IMAI

Copyright.
Edward Bausch.
1885.

D 211
B 35
1885

C. M. C.

TO THE READER.

It may seem to some persons an act of presumption for a maker of microscopes and microscope accessories to enter the field of authorship and attempt to supplement the valuable labors which in recent years have made the use of the microscope an indispensable aid in the advancement of science.

To such, if any, I submit, that being a producer of microscopes and their accessories, I have had opportunity to become acquainted with the lack of general knowledge of the fundamental principles of the instrument and the best method of *technique*, even among owners of microscopes. Indeed, with so many complications, with almost unlimited powers and uses of the instrument, the beginner cannot fail to feel the need of a guide and adviser.

In order to accomplish the greatest good, I have started out in this little *Manual* with the supposition that the purchaser, or owner, is a beginner and absolutely ignorant of the microscope and everything which pertains to it, and therefore have attempted to convey,

step by step, in as simple language as I could command, information which will, I trust, lead to ease of manipulation and give both pleasure and profit to those for whom it was specially written.

With these, its purposes and hopes, I beg for my self-imposed labor a friendly reception.

EDWARD BAUSCH.

CONTENTS.

SIMPLE MICROSCOPES.

| | |
|--------------------------------------------------------------------------------------------------------------------------------|---|
| Purpose of the Microscope; Kinds of Microscopes; Magnifying Power; Using Magnifiers; Aberrations; Achromatism, - - - - - | 7 |
|--------------------------------------------------------------------------------------------------------------------------------|---|

THE COMPOUND MICROSCOPE.

| | |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|
| Description of Parts; Base or Foot; Pillar; Arm; Body; Nose-Piece; Society-Screw; Eye-Piece, or Ocular; Draw-Tube; Collar; Coarse-Adjustment; Milled Heads; Fine-Adjustment; Stage; Clips; Cen- tering Screws; Mirror; Mirror-Bar; Sub-Stage; Sub-Stage Bar; Diaphragm; Optical Axis; Object; Slide; Cover-Glass, - - - - - | 12 |
|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|

OBJECTIVES AND EYE-PIECES.

| | |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|
| Objectives-Classes; Systems; Angular Aperture; Achro- matism; Resolving or Defining Power; Flatness of Field; Penetration; Working-Distance; Magnifying Power; Selecting Objectives; Eye-Pieces—Huyghen- ian; Solid; Periscopic; Flatness of Field; Size of Field, - - - - - | 17 |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|----|

REQUISITES FOR WORK.

| | |
|--------------------------------------------------------------------------------------|----|
| Working Table; Room; Light; Position of Light; Which Eye to Use; Order, - - - - - | 40 |
|--------------------------------------------------------------------------------------|----|

HOW TO WORK.

| | |
|-------------------------------------------------------------------------------------------------------|--|
| To Set Up the Instrument; Centering Stage; Illumina- tion; Attaching High Power Objectives; Double | |
|-------------------------------------------------------------------------------------------------------|--|

| | | | | | |
|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---|---|---|---|----|
| Nose Piece; How to Work; Dark Ground Illumination; Polarized Light; Cover-Glass; To Draw; Camera Lucida—How to Use It; Determining the Magnifying Power, | - | - | - | - | 45 |
| ADVANCED MANIPULATION. | | | | | |
| Dry-Adjustable Objectives; Immersion-Adjustable Objectives; Test-Plate; Immersion Objectives on Test-Plate, | - | - | - | - | 66 |
| SUB-STAGE ILLUMINATION. | | | | | |
| Objectives and Eye-Pieces; Hemispherical-Lens; Wenham Button; Woodward Prism; Narrow Angle Condensers; Wide Angle Condenser; Ward's General Illuminator; Large Sub-Stage Condenser, | - | - | - | - | 79 |
| CARE OF A MICROSCOPE. | | | | | |
| To Take Care of a Stand; To Take Care of Objectives and Eye-Pieces, | - | - | - | - | 85 |
| CONSIDERATIONS IN TESTING OBJECTIVES, | - | - | - | - | 89 |

SIMPLE MICROSCOPES.

Purpose of the Microscope.—The Microscope is an instrument which magnifies objects, so that we are better able to examine their structure than is possible with unassisted vision.

Kinds of Microscopes.—Microscopes may be divided into two classes—simple and compound—the difference between the two being that with the former, the object is viewed directly, while with the latter a magnified image is observed: while the first shows the objects in their true position, the latter shows them reversed, so that what is right in the object is left in the image, and when an object appears to be moving in a certain direction, the movement is in reality the reverse, and must be moved accordingly to keep it in view.

Magnifiers.—Simple microscopes are usually termed magnifiers, and, whether consisting of one or more lenses,



Fig. 1.

always remain simple. The most common are those with one or several double-convex lenses, (Fig. 1.)


The shorter the radii (the more curved the surfaces) are in these, the greater will be the magnifying power, and the higher this is, the less of the object's surface can be

seen at once. Each additional lens increases the magnifying power in proportion to its curvature. The distance between the lens and the object, when this is seen most distinctly, is called the *focus*; at the point where the object is most distinct, the lens is said to be *in focus*; when indistinct or blurred, *out of focus*.

Magnifying Power.—Unless a microscope is known to come from the hands of a reliable firm, any claim as to magnifying power should be accepted with reserve. In former years, when the country was over-run with cheap foreign productions, the most fanciful claims were made in this direction. It is evident that a lens magnifies an object equally in all directions; this is said to be in *areas*, and is the square of the *linear*, so that if an object is magnified 4 times in the linear, it is 16 times in area. The commonly accepted term to express magnifying power of simple, as well as compound microscopes, is in *diameters (linear)*. A single lens of 1 inch focus magnifies about ten diameters; one of 2 inch focus, about 5 diameters; one of $\frac{1}{2}$ inch focus, 20 diameters, and so on. In a lens of high magnifying power, the focus is ordinarily about twice the diameter, so that if a lens is $\frac{1}{2}$ inch diameter its focus is about 1 inch. This may, however, be more accurately determined by projecting, say a flame or window-frame, upon a white piece of paper; the distance between the paper to the center of the lens, when the image is most distinct, is its *focal distance*. When a lens is two inches or more in diameter, it is usually termed a *reading glass*.

Using Magnifiers.—In using magnifiers, the lens should be held close to the eye and such a position taken that the object will receive the best illumination. In the lenses of equally convex surfaces, it is immaterial which side is held toward the eye; but when plano-convex lenses are used, the plane side should always be toward the eye, as it gives the flattest *field*.

Aberrations.—Two factors arise which prevent the advantageous use of more than about 25 diameters in magnifiers; they are called the chromatic and spherical aberrations. The first is the term employed when the object is apparently fringed with color, predominantly blue and yellow; the second, when all but the central portion of the lens shows the object indistinctly; these faults increase with the magnifying power. In the case of a combination of several lenses, they may be partially overcome by interposing an opaque plate with a small opening, called a *diaphragm*, between them, which cuts off the outer or marginal rays, or the lenses may be made of a smaller diameter. An incision may also be cut into the glass equally between the two surfaces, when,

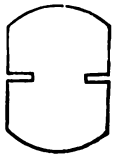

from the name of the inventor, it is called a *Coddington*. Fig. 2 shows a section of a Coddington, while Fig. 3 shows it in its mounting.



Fig. 3.

Achromatism.—The most approved method, however, for eliminating these appearances is by the use of one or two concave *flint glass* lenses in connection with the double convex *crown glass* lens. When the color or chromatic aberration is thus removed, the lenses are said to be *achromatic*, and when both the chromatic and spherical aberrations are avoided, the lens is called *aplanatic*, and is then said to be *corrected*. An achromatic lens, composed of one flint and one crown



Fig. 4.

glass lens, is called a *doublet*; (Fig. 4) one with two flint glass lenses and one of crown glass is called a *triplet*. Fig. 5.)

The latter is the best form, as it gives the highest correction; such a lens (it is thus called from the fact that the lenses are cemented together and act like one) may be held with either side toward the object with equally good results, and may also be held at quite an obliquity, without loss of definition; this feature is important, as it is almost impossible to give a lens a theoretically correct position to both the eye and object with the unaided hand.

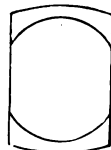


Fig. 5.

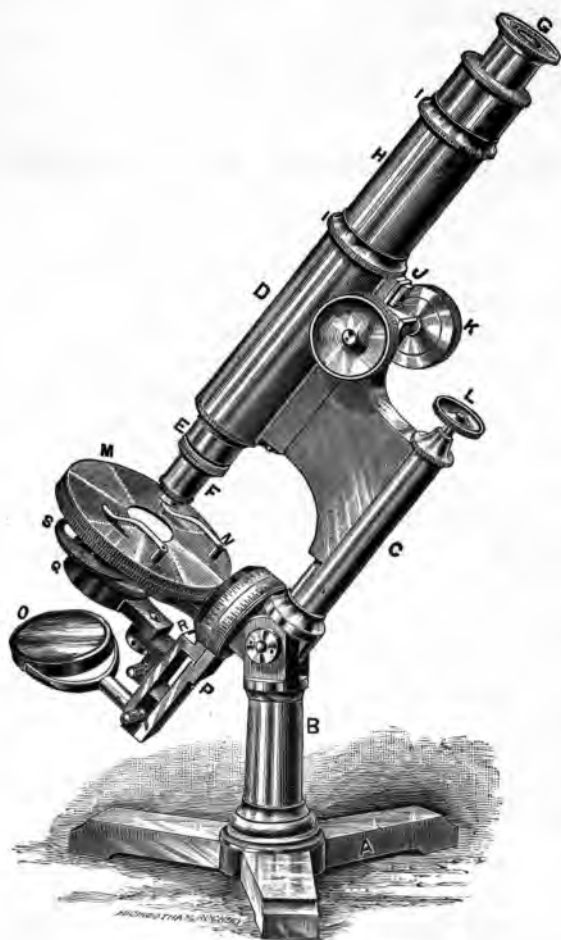


Fig. 6.

THE COMPOUND MICROSCOPE.

As was previously stated a magnified image is observed in the Compound Microscope. Any two lenses, one of short, the other of a long focus, placed sufficiently far apart, will attain this object, and this was for years the method of its construction. It was found that on account of the extreme magnifying power, mechanical appliances were required to keep the lenses steady during observations and at their proper distances, also that provision was necessary for adjustment and illumination. These adjuncts were so necessary that the instrument could not be conceived without them, and the entire apparatus is thus, by the force of usage, called a microscope, and the instrument without the optical parts, a *stand*.

Description of Parts.—As it is necessary for the student to become conversant with the terms of the various parts and to understand their use, we give an illustration (Fig. 6) with letters, and append a list giving the names.

A. Base or Foot.—This is the foundation of the instrument. It usually rests upon three points (or should do so) and is of such a weight that it keeps the instrument firm when it is in an upright

or inclined position. The revolving plate, when this is provided, by means of which the upper portion of the instrument is revolved, without changing the position of the base, is considered a part of it.

B. Pillar.—It is that portion which is fastened to the base and may be one or two, according to the construction of the stand. It carries upon its upper end the *joint* or *axis*.

C. Arm.—This is connected with the pillar by the joint and supports all the working parts of the instrument.

D. Body.—This is the tube-portion to which the optical parts are attached.

E. Nose-piece.—This is an extra piece which is attached to the lower part of the tube.

Society Screw.—This is a standard screw which is cut into the nose-piece, and is called so from the fact that it was first established by the Royal Microscopical Society of London. It is also called the *universal* screw, and is in general use in this country and England; it has lately been adopted by some firms on the Continent of Europe.

F. Objective.—This is screwed into the nose-piece and is called so because it is nearest the *object*. It is the most important of the two optical parts (of the microscope proper) and upon its perfection the distinctness of the image and therefore the value of the instrument almost entirely depends.

G. Eye-piece or Ocular.—It is called so because it is nearest the eye and is the remaining optical part. It magnifies the image given by the objective. This and objective will be treated more fully later on.

H. Draw-tube.—This is that portion of the body which moves in the outer sheath and which receives the eye-piece. It is provided for the purpose of attaining variations in magnifying power and as a matter of convenience while working.

I. Collar.—This is a ring which is attached to the draw-tube and is usually provided with a *milled edge*.

J. Coarse Adjustment.—This is a provision for moving the body quickly back and forth for adjusting the focus approximately. It is done by a sliding *rack* and stationary *pinion* (not shown in cut) or a sliding body in an outer sheath.

K. Milled Heads. These are attached to the shank of the pinion, which is revolved by means of them and are usually large to give sensitiveness to the movement.

L. Fine Adjustment.—This is slow moving and serves to get an exact focus. It is attained by a fine thread, provided with a milled head, and acts upon the body, either directly or by levers. This as well as the coarse adjustment should be extremely sensitive and should not have the least side or lateral motion. The fact that either of them have it, is evidence of poor workmanship.

M. Stage.—This is the portion on which the object is placed for examination. It is attached to the arm and may be either permanently fixed, or revolving or mechanical, the latter by moving the object by mechanical contrivances instead of by fingers. Metal and glass are used in its construction.

N. Clips.—These are two springs which are attached to the upper surface of the stage and serve to hold down the object.

Centering Screws.—These are provided for moving the stage in different directions to bring the center of its revolving motion in the center of the field.

O. Mirror.—This is used for reflecting and condensing light upon the object. As a rule two are used, one plain and the other concave. The first gives a comparatively weak light, while the second concentrates it and gives it more intensity.

P. Mirror-bar.—This carries the mirror and by a sliding arrangement allows the variations in distance of the mirror to the stage; it also swings in a circle around the object in order to illuminate it from any direction.

Q. Sub-stage.—This is a ring below the stage to receive various accessories which may be required. It is sometimes fixed to the stage but in the best instruments it is separated from it and is provided with an adjustment to vary its distance from the object.

R. Sub-stage bar.—This receives the sub-stage and permits its adjustment. This as well as the mirror-bar should be on an axis in the plane of the stage, so that whatever position they may be in, relative to the object, the distance from this to the sub-stage or mirror does not vary, except when made to do so.

S. Diaphragm.—This is a perforated, revolving disk, attached either to the stage or sub-stage. It has holes of different sizes so that the amount of light from the mirror may be modified.

Optical Axis.—This is an imaginary line which passes from the center of the eye-piece through the body, objective, stage and sub-stage to the mirror. Whatever lies in it is said to be *centered*.

Object.—Is that which is examined and placed upon a

Slide.—This is a thin plate of glass, generally 3 inches long by 1 inch wide.

Cover Glass.—This is an extremely thin piece of glass, round or square, which is placed upon the object, either for flattening or preserving it, or both.

OBJECTIVES AND EYE-PIECES.

Although considerable magnifying power may be attained by the use of two single lenses arranged in a compound form, there is no advantage in it, from the fact that the faults in the lenses are correspondingly magnified, and these are so considerable that they destroy what it is the purpose of the microscope to give—a distinct image.

Objectives, Classes.—Objectives may be divided into two classes, *dry* and *immersion*; in the former no intervening medium except air exists between the cover and objective, while in the latter a fluid is used to connect the upper surface of the cover to the front surface of the objective. The use of immersion fluid has several advantages, the first of which is that the objective may be made to give better performance, as will be explained later on; the second is that more light will be transmitted, as there is less loss of it by refraction. It must be observed that an objective which is made to use dry can not be used with immersion, nor *vice versa*.

There are two fluids in general use at the present time, *water* and *homogeneous fluid*. The latter expression means *of the same kind*, and refers to the fact that


the fluid has about the same refractive and dispersive power as glass, so that when this fluid fills up the space between the two surfaces of glass, a ray of light passes through the three mediums as if they were one body.

Objectives are sometimes called *powers*, and in this sense are divided into three classes: *low*, *medium* and *high*. Dr. Carpenter classifies them as follows: *low powers*, 3 inch, 2 inch, $1\frac{1}{2}$ inch, 1 inch, $\frac{2}{3}$ or $\frac{3}{4}$ inch; *medium powers*, $\frac{1}{10}$ inch, $\frac{1}{8}$ inch, $\frac{1}{4}$ inch, $\frac{1}{6}$ inch; *high powers*, $\frac{1}{8}$ inch, $\frac{1}{6}$ inch, $\frac{1}{10}$ inch, $\frac{1}{12}$ inch, $\frac{1}{16}$ inch, $\frac{1}{20}$ inch, $\frac{1}{24}$ inch.

As the objective is the most important of the two optical parts, it follows that this must be as free from faults as possible and all that human ingenuity and skill can devise is utilized to attain this end. The advance in the perfection of the objective has been step by step and each era, was at the time, considered by many authorities the limit to further improvement. Each advance was signalized by a marked opposition and disbelief of its possibility. It is therefore of inestimable credit to the pioneer objective-makers, and notably among these two Americans, who by quiet but stubborn application, disproved previous claims and opened the way to further improvements. A theoretical limit has been fixed on the capacity of the microscope, which, according to our present knowledge can not even be reached, but it is still safe to say that the end is not yet. A modern objective of the highest capacity may be considered a work

of art, and there are few productions of the human hand which exact so much untiring application, ingenuity and skill.

Systems.—An objective is said to consist of systems which may vary in number from one to four and five; two and three are however mainly in use. They are the individual portions consisting of one, two or three lenses, which when more than one, are

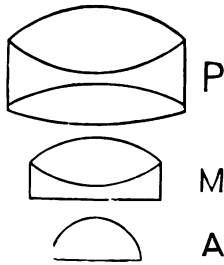

cemented together and make up the objective. An achromatic single system may consist of two or three lenses, and a three or four system objective may consist of as many as seven or eight lenses. The systems are called in their order: *anterior* or *front*, *middle* and *posterior*. When one

consists of two lenses it is called a *doublet*, when of three lenses a *triplet*. Thus in Fig. 7, A is the anterior, M the middle and P the posterior systems; thus also A is a single system, M a double and P a triple one.

The various features which must be considered as determining the quality of an objective are: angular aperture, achromatism, resolving or defining power, flatness of field, penetration, working distance and magnifying power. Although these attributes may be considered separately, some of them go hand in hand. The presence or extent of one necessarily involves or precludes another.

Angular Aperture.—The angle which the most extreme rays, which are transmitted through the objective, make at the point of focus, is called its angular aperture, or in short its *angle*, and of all the qualities in an ideal objective, this is the most important. Thus in Fig. 8, D is considered the point of focus, and C D E the angular aperture.

The above definition has its limitations, however. While in objectives of proper construction it holds true, there are many in which it is not the case. For instance, an objective may be so constructed that it may transmit a considerable number of rays in excess of those which combine to form an image, and it is evident that these should not be considered as belonging to them.

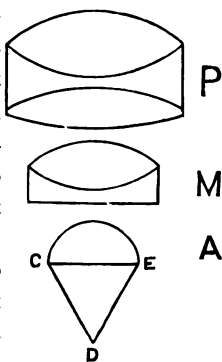


Fig. 8

As there are many objectives of the same power, but of different angular aperture, there are again others of varying power, but of the same angle. Other things being equal, it is the angular aperture of an objective which determines the quality. It is expressed in *degrees*, and is also spoken of as being wide, medium or narrow, although this is indefinite and depends considerably upon the power of the objective; while the angle may be excessively wide for a low power, it may be narrow for a higher one.

For many years the extent to which angular aperture could be carried was a matter of controversy, as was also the use of objectives of wide and narrow angles for different directions of work. It is, however, a matter of congratulation that the question is at rest, although it has served a good purpose in promulgating a better knowledge of the subject.

All objects emit rays, and it is evident that those coming from one point and contained in a large angle are more numerous than those in a small angle; also, that as the angle more nearly approaches 180° , the rays will not only be larger in number, but will embrace more of the surface of the point; if we imagine the point to be a sphere or globule, it is apparent that we see its upper half more distinctly, when in addition to a top view we also get a side view, than when we get a top

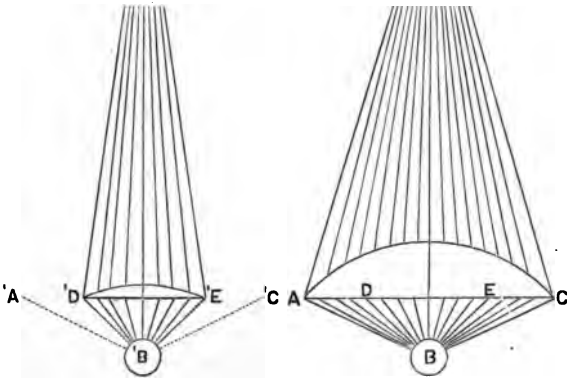


Fig. 9.

view only. As this is true with the unassisted eye, it is also a fact with the microscope. This is shown in Fig. 9. It is assumed that two objects, B and B', are equally bright, and therefore emit the same number of rays; for the purpose in hand, it is sufficient to consider only those which reach the plane surface of the large lens; if an equal space is imagined over B', the same number of rays will be contained in this; therefore, the cone contained in the angle A' B' C' will contain as many rays as that contained in A B C; but as the lens D' E' is considerably smaller than A C, only as many rays can enter it as are contained in the angle D' B' E'. As the rays contained in the angle A B C and D' B' E' are carried through the two lenses, which are supposed to be of the same magnifying power, the image formed by A C will be considerably brighter than that formed by D' E', and will therefore show more of its structure, as will be shown hereafter. Besides this, the rays comprised in the angle A B C cover more of the object's surface than do those of the angle D' B' E'.

It is therefore indisputable that objectives of the same angular aperture, but of different magnifying power (within ordinary limits), will show the object equally well, provided they are otherwise of the same quality, and also true that in objectives of the same power but unequal angular aperture, the one of wider angle will show an object more brilliantly than the other, and, if the difference be considerable, will show structure of which no trace can be found with the narrow angle. These are facts which are based upon natural laws, but there are

other conditions to be considered in connection with them, which will be treated hereafter.

It very often happens that objectives from different makers, but of the same angle, show a considerable variation ; this does not prove that the above principle is wrong, but is evidence that greater care or skill has been bestowed on one than on the other.

By arranging an objective with an *immersion front*, its angular aperture may be considerably increased over that of a *dry front*, and this explains why better results may be obtained with the former than with the latter.

Achromatism.—As has been stated before, when single lenses are made to give a high magnifying power, the chromatic and spherical aberrations prevent corresponding advantages ; and as the objective gives the image which is magnified by the eye-piece, it is evident that if they exist in it, they are increased by the ocular, and that especial care must be given to exclude all faults as much as possible from it. Even with the use of flint glass, it is impossible to free the objective entirely from color ; there will remain a residue of green and purple, and these colors will fringe the object. These are called the *secondary spectrum*, and their presence in an objective is usually evidence of the highest correction.

The amount of color in an objective depends somewhat upon the power of the eye-piece, and becomes more visible as a higher power is used. Color outside of the secondary spectrum is not always

prejudicial to an objective ; for, if in two, one shows the structure of an object with a slight amount of it, the other does not show the structure but gives a nearly colorless image, it goes without argument to say that the first gives the best results and is therefore preferable.

If on increasing the distance between the objective and object the latter shows a marked bluish color, and when the distance is decreased a yellowish-red color, the objective is chromatically *under-corrected* ; if, however, the conditions are reversed—if the object shows a yellow color when the distance between it and the objective is increased, and blue when decreased, it is *over-corrected*. If the colors are not sufficiently pronounced to judge the corrections, they may be determined more positively by means of oblique light. If, when the mirror is swung to the right, the right edge of the object is yellow and the left edge blue, the object is chromatically *over-corrected* ; if the right edge is blue and the left yellow, it is *under-corrected*.

Resolving or Defining Power.—This is the quality in an objective by which we are enabled to see the intricate structure and finer details in an object. It depends upon the correctness of the chromatic and spherical aberrations, upon angular aperture, and of course upon the perfection of the mechanical work. The power of resolving in an objective is indicative of the perfection of the micro-

scope, for it is almost entirely dependent upon it for its quality.

When an objective is said to *resolve* a structure or a certain number of lines, it means that it shows them under certain conditions of light. It may resolve *easily* or only *glimpse* them—the latter when they are hardly to be distinguished. Every degree added to the angular aperture of an objective increases its resolving power, and the theoretical capacity of every degree has been mathematically determined. However, this standard is only reached approximately and to a varying extent. It is not by any means said that every objective of a certain angular aperture will have a corresponding resolving power; it is at this point that the acute accuracy of work and superior judgment of proper corrections will invariably give the best results.

It is an error to suppose that the resolving power may be improved by merely increasing the magnifying power. It is an invariable quality of an objective and has a fixed limit. The extent to which it may be approached depends upon the nicety of manipulation, but no amount of increase in magnifying power by the eye-piece or any other means will carry it beyond it; on the contrary, it will lose in this respect if carried beyond a certain point.

Flatness of Field.—The *field* in a microscope is that portion which is observed in the eye piece, and its flatness may be observed when focused on a flat object—preferably a micrometer. It is said to

be flat when all portions of the object are seen over the entire field at once without further focusing. When not flat, it will be found that as the image approaches the edge of the field it becomes more and more indistinct, and that the objective must be

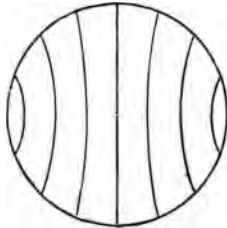


Fig. 10.

correspondingly adjusted; in many cases it remains indistinct or blurred, and this may be considered the most serious fault. In the case of looking at straight parallel lines, such as in a micrometer, they will appear to become more curved as they near the edge, as shown in Fig. 10.

Flatness of field mainly depends upon the correction of the spherical aberrations, and as under the best conditions the latter cannot be entirely eliminated, it is impossible to attain absolute flatness, except with eye-pieces especially made for this purpose. It may, however, also be due to a faulty eye-piece; in this case it can fairly be determined, by observing whether it shows equally in different objectives. With beginners, especially, it is usually most complained of, owing probably to the fact that it is the most easily noticeable. Although a desirable quality, we consider it a matter of minor importance, and in the choice between two objectives, in one of which the predominant feature is resolving power and in the other flatness of field, our choice would invariably be for the former.

Penetration.—This is the quality which enables us to look into an object—to observe different planes at one time. In the mind of the writer, it is of no special importance, or, at any rate, not as much as is claimed for it, and, if desired, is easily attained. It depends upon magnifying power and angular aperture, and decreases with the increase of either of these. Objectives are generally not constructed with any reference to it; it is a natural consequence of certain conditions.

Penetration and resolving power are antagonistic, or, at any rate, in an inverse ratio, and can only be combined to a certain extent. In two objectives of the same power and angle, one cannot have penetration as a special feature and the other resolving power; they will be almost similar in these qualities, provided that they are similarly corrected. However, if they are not similar in their angular aperture, the one of narrow angle will have more penetration than the other. In objectives of the same angle but different power, the one of low power will have in itself more penetration; it will be similar in its action to the eye, which, when an object is close to it, can distinguish but one portion of it distinctly, while, as its distance to the eye is increased, can distinguish various parts of it lying at different distances, and will finally see other objects outside of it. By looking at an object at 5 feet distance, only this can be seen plainly; but, at 10 feet, others quite a distance in front or back of it can be seen as well.

Working-Distance.—This term, strictly considered, is an invariable quality of the objective, and is the distance between the front lens in the objective and an uncovered object, when the objective is in focus and is corrected for that object. All objectives require a certain amount of projecting metal to protect the front lens, and this with a certain thickness of the cover-glass lessens it. In objectives with fixed mountings this may be, and with thick cover-glasses is considerable. As it is comparatively unimportant, however, for the working microscopist to know the working distance *per se* of his objectives, but of considerable moment to know what the actual space between the objective and cover-glass is, it would be well, in the mind of the writer, to express it as *available* working distance.

In objectives of low and medium power, it is of little consideration; but where it must be expressed in 1-100 or 1-1,000 inch, it becomes a matter of importance.

Working-distance is spoken of as being *long* or *short*, and varies not so much with the power as with the angular aperture; generally the working-distance decreases with the increase in angular aperture, and becomes greater as the angle becomes narrower; it was for a long time considered that these two properties varied according to a fixed rule, but this at the present time is not considered to be the case. While in objectives of the same angular aperture it may vary considerably, it may in others

of different angles be so that the one of wider angle may have the greatest working distance. The skill of the optician must in a considerable manner determine the amount of it.

It will be seen from the above that working-distance stands in no direct relation to the focal distance of the objective, neither to its nomenclature or rating, and, it may be added, that it is never as great as the focal distance of a single lens of the same magnifying power.

As may be imagined, there are a variety of opinions as to what constitutes long or short working-distance in a certain objective. No definite rule can be laid down for this, as it is conditioned by the skill and requirements of the manipulator. Although it is an important factor, the idea that it should in all cases be as great as possible, is erroneous, for, while it may be true in a dry objective, it may be the cause of annoyance in one with immersion. On several occasions it occurred in the experience of the writer that after an objective had been completed, it was found that its working distance was so large that the immersion fluid would run out from between the objective and cover-glass when the instrument was inclined, and it was necessary to change the objective with a view to decreasing its working distance, in order to allow its convenient use.

Magnifying Power.—This is a question of vital importance in a microscope, not so much as a quality for itself, as in connection with the resolving

power. The inquiry should not be simply how many diameters an instrument will magnify, but what the precision and extent of its definitions are under a certain magnifying power. If a high magnifying power is all that is desired, this may be obtained to an almost unlimited extent by means of simple lenses which may be procured at a small pecuniary outlay; but these do not give a distinct image nor do they make structure visible, which, be it remembered, is the purpose of the microscope to do.

The normal eye can distinguish from 200 to 250 lines to the inch, and in a microscope such magnifying power should be used, which will apparently bring the structure which is sought after *at least* up to this figure. In illustration, take a $\frac{1}{2}$ inch objective of 98° and a $1\frac{1}{2}$ inch eye-piece. An objective of this kind, properly corrected, resolves *pleurosigma angulatum*, in which the average lines are 60,000 to the inch. With the above eye-piece it is utterly impossible to see them, while if it is re-placed by a $\frac{3}{4}$ inch or $\frac{1}{2}$ inch, they can easily be distinguished. This is not owing to any peculiar quality of the eye-piece, but merely to the fact that by increasing the magnifying power, the dimensions of the object have been increased to such an extent that the lines have apparently been separated and become visible to the eye.

Beginners as a rule are apt to use too much magnifying power or *amplification*, and often attempt

to view a large surface with an objective which will show but a small part of it. It must not be forgotten that the apparent field of view is decreased as higher powers are used, and that a low power will give a better impression of a large, coarse object and its relative parts, not only because it makes a larger surface visible, but because it has more penetration.

In objectives of the same power, but of different angular aperture, the magnifying power and field will always be the same.

The following table, which has been compiled, will probably be of assistance to the beginner. After he has become better acquainted with his instrument his judgment will dictate to him what to do.

A power of 25 diameters will show a surface of about 1-5 inch diameter.

A power of 50 diameters will show a surface of about 1-10 inch diameter.

A power of 100 diameters will show a surface of about 1-20 inch diameter.

A power of 500 diameters will show a surface of about 1-100 inch diameter.

A power of 1,000 diameters will show a surface of about 1-200 inch diameter.

This table is approximately correct with a Huyghenian eye piece; with a Periscopic almost double the amount of such surface will be shown.

Magnifying power may be obtained by the eye-piece or objective, and the desirability of using one or the other for this purpose was for many years a matter of spirited discussion.

Objectives of the same angular aperture, but of different power, will give identical results, by bringing them up to the same magnifying power, unless the difference is considerable. In both objectives and eye-pieces the lenses decrease in size with the increase in power and consequently give less light; and while this one objection exists in the objective an additional one occurs in the eye-piece, in that the eye must be brought closer to the eye-lens and must be kept more strictly in the optical axis, which at a long sitting becomes fatiguing. By almost common consent a higher power than $\frac{1}{2}$ inch eye-piece should not be used in prolonged work, and between this and the lowest it becomes a matter of individual choice. All responsible manufacturers and dealers make up such outfits of stands, objectives and eye-pieces, which experience has shown are most generally useful.

It is a safe rule to follow in all work on recognized forms (objects of which the structure is known) not to use a higher power than is necessary to properly study them.

Selecting Objectives.—It is supposed that what has been stated regarding the properties of an objective, has been sufficiently explicit to enable the

beginner to make a suitable selection, when he finds that other objectives than those which he possesses are necessary. Although from a theoretical standpoint it is true that objectives are optically more perfect as they increase in their angular apertures, it is not by any means said that this is a rule in practice. A great portion of every-day work can be done with less loss of time, with more comfort and with the attainment of as high results with objectives of medium angular aperture and magnifying power. A 1.5 inch of 75° or 110° for instance is a very valuable objective in any outfit, although its performance is limited; when more is desired than this will give, recourse must be had, first to objectives with wider angular aperture (this should be considerable); second to those with more magnifying power.

In this connection it is considered important to mention the suggestions of Dr. Geo. E. Blackham, who from his knowledge of the optical principles involved in the microscope, and from his long experience and delicacy in manipulation, is pre-eminently fitted to make recommendations.

His ideal series is as follows:

OBJECTIVES.

One 4 inch 12° air angle = 0.10 numerical aperture.

One 1 inch 30° air angle = 0.26 numerical aperture.

One 1.6 inch 140° air angle = 0.94 numerical aperture.

One 1.8 homogeneous immersion 1.42 numerical aperture.

EYE-PIECES.

One 2 inch Huyghenian.

One 1 " "

One $\frac{3}{4}$ " "

One $\frac{1}{2}$ " Solid.

The 4 inch and 1 inch objectives are of course non-adjustable, while the two latter have cover correction.

This series is supposed to cover every direction of investigation. In a great many cases, of course, such a low power as would be given by the 4 inch, would be useless; in the above table the highest power is obtained by the use of a $\frac{1}{8}$ inch objective and $\frac{1}{2}$ inch eye-piece, whereas the same can be accomplished by a 1.12 inch objective and $\frac{3}{4}$ inch eye-piece, or 1.16 inch objective and 1 inch eye-piece. As stated before, this becomes a matter of individual choice.

Be distrustful of all objectives in which the chief inducement is a very low price; they are almost invariably useless for any reliable work. Have a distrust especially of all "nameless" objectives. It is safe to assume that if the maker can not attach

his name, he is doubtful of their superiority, if not convinced of their inferiority. Any maker of responsibility will say without hesitation, that he can produce objectives at less than one-half their present cost, if he had the assurance that they would be accepted as first put together, as the cost of merely making and mounting lenses is considerable less than the cost of making proper corrections. In this case however they would be of varying and inferior quality.

In purchasing a microscope a beginner may be easily misled by the enticing appearance of an object, which may be due not so much to the instrument as to the object itself, and if the optical parts are inferior, it will require but a short experience to become convinced of it—usually as soon as a comparison can be made with reliable work. The investment in one of these objectives is not only a source of disappointment, but usually proves to be a pecuniary loss, as it is usually followed by a fresh outlay in responsible work.

It is of ordinary occurrence that such objectives as have just been spoken of are sent to the writer's firm with the request to examine them and rectify the faults; but an examination almost invariably proves that the cost of doing this is considerably greater than purchasing a new objective of the same power, and it would not then even be equal to the latter.

Eye-Piece—Huyghenian.—This is now in general use, and consists of two plane-convex lenses.

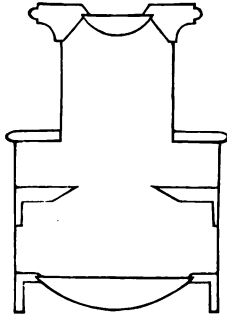


Fig. 11.

It receives its name from the inventor, who first applied it to the telescope. The *eye-lens* is the small lens nearest the eye, and the *field-lens* or *collective*, as it is also called, is the large one nearest the objective. A *diaphragm* is placed between them, and gives a sharply defined *field*. This eye-piece is also called *negative*, as its focal

point is between the two lenses (at the diaphragm) in contradistinction to a *positive*, in which the focal point is outside of and below the field lens.

Solid Eye-Piece.—This was the invention of the late R. B. Tolles, and also belongs to the class of negative eye-pieces. It is called *solid* from the fact that, instead of being composed of two lenses, it consists of one piece of glass, which is cut to a cylindrical form, and on the ends of which the proper curvatures are ground; the diaphragm is made by cutting a circular groove into the glass at the proper distance between the two surfaces, which is then filled up with an opaque pigment.

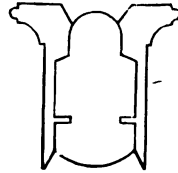


Fig. 12.

These eye-pieces are only made in high powers, as optical glass is usually not of sufficient homogeneity to make low powers, and their cost would be too considerable, without a corresponding advantage. For high powers they are superior to the Huyghenian, in that they give a better illuminated field, as there is less loss of light by absorption through the glass than by refraction at the two additional surfaces of the eye-lens and field-lens in the Huyghenian.

Periscopic Eye-Piece.—This consists of a triple eye-lens and single field-lens. Its predominant feature

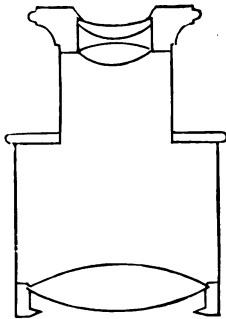


Fig. 13.

is a very large and flat field, with almost all objectives. In this respect it has a considerable advantage over the Huyghenian and Solid. It is positive and therefore well adapted for micrometer work, as it is focused like a magnifier, and its magnifying power remains constant, while with the Huyghenian it is variable, from the fact that the eye-lens alone is focused, thus varying its distance from the field-lens, and consequently the magnifying power.

Nomenclature.—The rating of eye pieces was formerly, and is to a considerable extent to-day, by letters. This method, however, is arbitrary, as the letters of different makers have a totally different

significance, so that nothing like a standard exists. This fact induced the American Society of Microscopists to endeavor to establish a universal method, and after the matter had been given careful attention for several successive years, it finally adopted the method, which rates them according to their magnifying power, the same as that which has been used in objectives. This gives an approximate idea of the magnifying powers; thus, an eye piece marked 1 inch or by a letter signifying the same, shows that it magnifies about 10 diameters; one of $\frac{1}{2}$ inch, 20 diameters, and so on.

Flatness of Field.—Although this depends mainly upon the objective, the absence of it may be owing to a faulty construction of the eye-piece. If it is so prominent as to be easily noticeable, and to the same degree with a number of objectives, it may be ascribed to the eye-piece. It must, however, be remembered that an absolutely flat field has not yet been obtained; it may be closely approached by decreasing the diameter of field to less than its normal size.

Size of Field.—Quite a general but erroneous idea prevails that the size of the tube has an influence on the size of the field. Except in eye-pieces of very low power, or with tubes of smaller than usual dimensions, this is not so. It must be remembered that a Huyghenian eye-piece admits of a definite size of field, and this is regulated by the opening in

the diaphragm ; the same size of opening is used in all of the same power, whether it is for an eye-piece of a large or a small diameter.

A misconception also exists as to the definition of field. Such inquiries are often made as: "as we understand it, a wide-angle objective gives a larger field"; but it does nothing of the kind. The angular aperture has no bearing whatever on the size of the field. The *field of view*, or that which is shown of the object's surface, is determined by the power of the objective and eye-piece.

REQUISITES FOR WORK.

The necessary material for preparing and mounting specimens will be found enumerated in books which are devoted to this purpose. It is the intention to make such recommendations in this chapter which, if not absolutely necessary, will be found convenient and will aid in facilitating work.

One of the first requisites for the proper use of the microscope is a thorough knowledge of its parts and an acquaintance with the optical principles involved. For this purpose the writer earnestly requests a perusal of the preceding pages, and is convinced that in cases where no previous knowledge of the instrument has existed, work will be done with far more ease, in much less time, and with a greater degree of satisfaction. Ignorance of the instrument's capacity may lead to an idea that it is inferior and thus be the means of its final abandonment; and in place of the anticipated pleasure there may arise a feeling of bitterness and disappointment for all future with everything connected with it. There are innumerable cases of this kind and they have induced a belief that it is difficult to acquire a practical manipulation of the microscope, whereas such is not the case when a limited time, properly applied, is devoted to it.

Working Table.—A firm table should be used, preferably one with three legs, as this will always be firm, no matter how uneven the floor is, and if it can be arranged, should be devoted to this purpose only. One with a round or square top of three feet provides ample room. Although not necessary, a table with a revolving top, provided with clamp, is very convenient, as with this two or more persons may make observations without changing their seats.

A very neat arrangement for a table-top is that suggested and used by Dr. J. E. Reeves. He places upon an ordinary table three or four thicknesses of white paper and upon these, a plate of polished glass, as large as the top; this can be procured of almost any glazier at a low price. It is pleasant to work upon and will not soil.

As in almost all cities there is more or less continual vibration from wagons upon the paved streets, the writer suggests an effectual remedy. Take a thin board, say half an inch thick, of a sufficient size to receive the microscope; fasten on the upper side at two opposite ends, cleats of about 1 inch square and countersink into these through the board four spiral springs of such tension that when they bear the weight of the instrument, the bottom of the board will be about $\frac{1}{4}$ inch from the table.

Have the working-table provided with drawers and arrange receptacles for the accessories, secure from dust, but at a convenient point to reach. When

the instrument is not in use, put it into its case or cover it in a manner, so that it shall be free from dust. For this purpose a large bell glass is best.

Room—If possible a room should be selected facing the north, as the light in this direction is most constant. It will prove a great saving of time if all or a portion of it can be permanently arranged to receive the entire working outfit. It should also be chosen with a view to its being free from disturbance.

Light. As stated, the light from the northern sky is most desirable, and that from a white cloud is preferable to that from a blue sky. On account of its intensity, direct sunlight should never be used; but if modified by a white curtain or reflected from a white wall it is excellent.

For lamp light an ordinary flat wick kerosene or student lamp is well adapted. The Hitchcock lamp, from its better combustion is still better, as its color more nearly approaches white. The ideal artificial light is that from an electric light. Gas light is not desirable as it is seldom sufficiently steady.

Position of Light.—The relation of the microscope to the source of light is principally a matter of personal convenience. With daylight it makes little difference whether it is at the front or side of the instrument, although the writer prefers it at the front, as the manipulation of the object does

not obstruct it; but the lamp should be placed at the right or left side, within easy reach of the hand for the purpose of controlling it. The writer suggests that the beginner make it a habit at the outset to place it on the side of the instrument opposite to the unoccupied eye, as the tube then places the latter in the shadow.

Which Eye to Use.—In a binocular instrument both eyes are used, but in a monocular only one is used, and it depends upon a trial which is best suited. A large proportion of persons are afflicted with astigmatism, often without knowing it, and when this exists it may be in one eye or when in both, may be to a greater extent in one than in the other. Its presence may prevent the eye from observing fine detail; but whichever eye is found to be best suited should be used. When both eyes are alike, it is sometimes advisable to change from one to the other.

It should be made a habit at the outset, and strictly adhered to, to *keep both eyes open*. A little difficulty may be found to do this, as the eye which is free will probably observe the objects upon the table; but as soon as the mind becomes fixed upon what it sees in the microscope, this impression disappears. After a time it will be found to require no exertion and will certainly add to the ease and comfort of the manipulator while working.

The Ward Eye Shade, Fig. 14, will prove of assistance in acquiring the above mentioned habit, and besides this, effectually excludes the light from the eyes. It is made of hard rubber and is attached to the tube of the microscope.



Order.—Among the requisites for successfully prosecuting work with the microscope are a strict observance of the instructions, even if they appear superfluous, a systematic way of doing work, and cleanliness. Have a place for every article which is required, so that the hand may immediately be placed upon it; after it has been used, clean it before putting it aside; keep strange hands from your apparatus, unless you are assured that a knowledge of its manipulation exists.

HOW TO WORK.

To Set up the instrument.—Draw the instrument from the case by grasping the base, and free it from dust. If it has *draw tube*, bring it to its standard length, which is indicated by a ring, by as little of a screw motion as possible; if the draw tube is highly polished, its surface will be best retained by observing the above precaution. See that the *mirror* and largest aperture of the *diaphragm* are in a central position—in the *optical axis*. After being convinced that the *eye-pieces* are clean, place one into the *tube*; then remove the *objectives* from their cases and after first having increased the distance between the *stage* and *tube* by means of the *milled heads* of the *pinion*, attach the lowest power to the *nose-piece*, by using both hands, being careful that it is as near as possible in the *optical axis* while screwing it on. Then incline the *body* by placing the left hand upon the *base* and drawing with the right hand upon the *arm*; be careful not to pull on the tube, as it may prove too heavy a strain upon this or the fittings. Incline the body of the microscope until the eye-piece about reaches the level of the eye, so that when an observation is made the position is as comfortable

as possible; the neck should not be strained, neither should the chest be compressed. Next place the *slide* with a transparent *object* upon the *stage*, by sliding it under the spring *clips* and get it as near as possible in the center of the opening; for an object anything near at hand, such as a piece of printed paper or cotton fibres will do. Watching the slide, adjust the mirror until it is seen that the light strikes the object; incline the head to the level of the stage, and observing the objective, rack it down to within $\frac{1}{4}$ inch of the object. Again placing the eye at the eye-piece, reverse the motion of the milled heads and observing the *field* continue the upward motion of the body until the image of the object appears in view.

Centering Stage.—If the microscope has a revolving stage, turn this to see whether the object or portion of it lying in the center of the field, remains in the optical axis. It was true when the instrument was shipped, but may have changed during transportation. If not centered, loosen the screw holding it to the arm, by means of the steel pin, just sufficiently that by the exertion of a little pressure it can be moved. After having observed first which portion of the object remains stationary during its revolution (this evidently is its center) move the stage so that this point will be in the center of the field, and then tighten the screw. If the point lies outside of the limits of the field, its direction can be noted and the stage moved accordingly.

Where *centering screws* are provided in the instrument, this is a simple matter.

If the *coarse adjustment* does not prove sensitive enough to *focus* easily, adjust by the fine adjustment by taking the head of the *micrometer screw* between the thumb and first finger and move toward the right or left as may be necessary.

It may here be said in passing that the rack and pinion should be so well fitted that they should permit the adjustment of low power objectives with the greatest ease, and should work *without the slightest lost motion* with a 1-5 or 1-8 inch objective. This point is the criterion of workmanship in an instrument, and if it is found to have the least *back-lash*, or is not perfectly smooth, it may safely be assumed that the instrument is of inferior workmanship.

If the fine adjustment does not act, the screw has either come to its stop or has "run out," and must be brought into action again; the range of movement in almost all fine adjustments is quite short and constant care must be taken to keep it at about a medium point. If the object is found not to give a full view or is not in the center of the field, it must be moved on the stage, but it must be remembered that a movement in one direction causes an apparent opposite movement in the field. At first this movement will be in jerks, but after a little practice the necessary sensitiveness of touch is acquired to give it more steadiness.



Fig. 22.

and is attached to the microscope in the same manner as an objective. In those made by the Bausch & Lomb Optical Co., their 3-4 inch and 1-5 inch, 1 inch and 1-4 inch objectives, student series are fitted in pairs so that they correspond in focus without adjusting and are exactly centered. Each end of the nose-piece is marked for the proper objective.

How to Work.—It is now supposed that the instrument is ready for work. To start, it is well for the beginner to provide a few prepared specimens, as these will help him considerably, if it is his intention, as it should be, to prepare them later himself. Whatever branch of study he is going to follow, a slide of *Pleurosigma angulatum*, dry, will be valuable to practice upon and to determine the quality of his higher power objective. In this latter respect, however, the writer would advise the beginner to guard against expressing an opinion too soon. He knows of many cases where the optician's claims were flatly denied, when often a few words

of advice by letter or a few minutes of manipulation would resolve the diatom, when it had not been done by the beginner, and would do it so easily, that it became a wonder how it could be avoided in any hands.

For a low-power objective, the *proboscis* of a *blow-fly* is probably the most suitable and at the same time, most interesting object. Place this upon the stage, and, after getting it as close as possible to the center of the opening in it, focus by means of the coarse adjustment. If only a portion of it can be seen and if it is desired to see a larger surface, the length of tube may be contracted by means of the draw-tube. In this case the object will be placed *out of focus*, and another adjustment becomes necessary. If a higher power is desired, the draw-tube may be extended. Observe whether the field is well illuminated, and if not, bear in mind what has been said on this subject. If the object appears milky or the light is so intense as to be painful to the eye, which is of usual occurrence to the beginner, the diaphragm should be turned from one aperture to another until a marked difference is seen; or, the plane mirror should be used. In this connection, it is well to state that the above precautions should always be observed with low powers, unless the object is thick. Now use the micrometer screw and note carefully the beautiful structure which is opened to view. After

sufficient time has been spent upon this, the objective may be replaced by a higher power and the object by a slide of *P. angulatum*; focus upon this, being mindful of the suggestions previously given, and do not fail to observe what has been said in regard to a well illuminated field. If lamp-light from a flat wick is used, turn the edge of the flame toward the mirror, and use the concave side of the latter. If the diaphragm is in an adjustable sub stage, bring it as close to the stage as possible, or, whether here or attached to the stage, it may as well be removed for the present. Observe now whether outside of the central rib any lines can be seen upon the surface of the diatoms; if not, vary the distance of the mirror from the object; or, if lamp-light is used, bring the lamp closer to or remove it from the instrument in one line, so that the illumination will not disappear. If this does not bring out the lines swing the mirror-bar from the *central* position into an *oblique* one, on the side opposite to that of the light and readjust the mirror; in doing this grasp the ends of the *mirror-fork* between the thumb and middle finger and move the mirror by the first finger. If the field can not be evenly illuminated, it is evidence that the mirror is beyond the limit of angular aperture in the objective, and it must therefore be brought back until it is. It must here also be noticed that if the diaphragm is still attached to the instrument and does not swing with the mirror, it may also be the

means of cutting off light. By means of the micrometer screw carry the fine adjustment back and forth beyond the plane of the object and observe closely whether any lines can be distinguished. It is very probable that they will show ; but if not, the cause should be determined. It may be that the magnifying power is not sufficiently great, and in this case a higher power eye-piece should be used, or the cover glass may be more or less than the normal thickness, which would cause a spherical over or under-correction in the objective. In this case the lines would appear when the diatom is not in focus. If the objective is a non-adjustable one, the proper correction may be approximately reached by means of the draw-tube. If the lines appear over the plane of the object, it shows over-correction, and the length of tube should then be decreased, and contrary when the lines show below or beyond the plane of the object. If the above directions have been followed, the lines cannot fail to be seen with a moderately good 1-4 or 1-5 inch objective ; but if they are not, the trial should be repeated. Again, be careful to have no obstruction between the course of rays from the mirror to the stage ; get good illumination on the object ; observe well ; keep the instrument in such a position that the object is not illuminated from any other direction than from the mirror.

When the diatoms are *resolved* in this manner, the lines will appear to be *diagonal* in some ; *longi-*

tudinal or *transverse* in others, according to their position; and, if the resolution is very good, these lines will further resolve themselves in minute beads of a hexagonal form.

It will now be well to bring the mirror more nearly to a central position; do this at intervals of about 10 degrees, and note the appearance at each decrease of obliquity. It will be found that as the mirror approaches the optical axis the lines will appear to become more faint, and may disappear before central illumination is reached; in this case it will be well to begin again. An endeavor should be made to make each attempt give better results than the preceding one. Repeated trials will not only impress the various phenomena upon the mind, but will cause a notable improvement in manipulative skill, and thus a better performance in the objective.

Until now we have assumed that *transmitted* light has been used. We will now suppose that the object is not sufficiently transparent to use this method; the object is then said to be *opaque*, and requires a different procedure. We will say that it is desired to examine an insect; it may be attached to a slide, or, what is better still, may be fastened in a stage forceps—(Fig. 17),—as it may then be turned and viewed from all sides. The low power objective should again be attached; after having been focused, it will be found that the light is insuf-



Fig. 17.

ficient to illuminate it. The mirror-bar should now be swung

upon its axis around the stage to a point above it, so it will be at an angle of about 45 degrees to its surface. If a lamp is used and is in the same position as when used with transmitted light, it is probable that the tube of the instrument will obstruct the light, and it is then well to move it toward the front. Using the concave mirror, adjust it so that the light will be concentrated upon the object, by watching it directly, and then observe through the tube. If it is not sufficiently illuminated, continue to adjust the mirror; also vary its distance from the object and swing the mirror-bar to a higher or lower point. It often occurs that, under the best conditions, the need of better illumination is felt; in this case a *bull's-eye condenser*—(Fig. 18)—should be procured. It will be found that this will become a useful and perhaps necessary accessory in work outside of this. Place it close to the instrument and set the bull's-eye between the object and source of light, with the *plane* side toward the object; if an ordinary hand lamp is used, it will be necessary to elevate this to about the height of the eye-piece, and if it is to be used often in this position, a special support should be made for this purpose.



Fig. 18.

Low power objectives are usually used on opaque objects, but sometimes a higher power is desired. Unless one is constructed with a view to opaque illumination its working distance is usually so short that it will prevent the light from striking the object. A 1-4 or 1-5 objective, of 75 degrees, has sufficient working-distance, and its mounting is made conical in the front, so that it will allow it.

Dark Ground Illumination.—This method is not in general use, probably because it requires a special accessory, although it yields beautiful effects. It is accomplished by means of a paraboloid—(Fig.

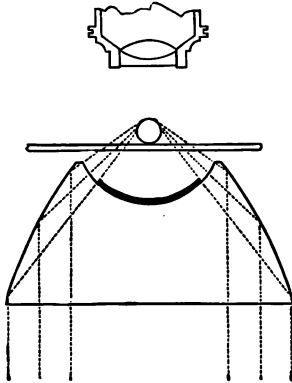


Fig. 19.

19)—which is attached to the sub-stage. As will be noticed in the illustration, the lower surface of the paraboloid is plane, and the light passes through this without undergoing any change. When it reaches the polished parabolic surface it is reflected to one point, according to the simple optical law that the angle of reflection is equal

to the angle of incidence. An opaque stop, which is cemented to the concave surface, prevents the light from passing through the central portion of the paraboloid. The object is thus illuminated on all sides by such an obliquity of light that it does not pass into the objective; the object stands out in relief, pleasantly illuminated on a dark back-ground. In using the paraboloid, the plane mirror should be used, and it is necessary to vary its distance from the object in order to attain the best results.

Polarized Light.—There is still another method of illumination. The light from the mirror is *polar-*

ized (undergoes a change) by passing through two *Nicol* prisms, one of which is fastened into a revolvable mounting in the sub-stage, and is called a *polarizer*. (Fig. 20.) The other is in a mounting above



Fig. 20.



Fig. 21.

the objective is called an *analyzer*;—(Fig. 21)—both in conjunction are called a *polariscope*. Although this is not an accessory for general purposes, it is nevertheless one which is very popular, and produces the most striking display of brilliant colors imaginable. Although it is strictly a scientific apparatus for use on crystals, it may be made the means of passing extremely pleasant and profitable hours. Any object in a crystalline form is suited, also sections of rock. Another interesting study with this apparatus is the observance of the crystallization of table-salt, as the water evaporates, after it was first dissolved.

Cover-Glass.—Thus far no attention has been given to the use of the cover-glass, although it is an important factor in reaching good results. In pre-

liminary examinations of solid objects with low powers it may be dispensed with; but where fluids are used, whether with low, medium, or high powers, it *should always be used*. A drop or small quantity of fluid placed upon a slide assumes a spherical form, and, on viewing it with a low power, it will be found to give a distorted field, and will cause disagreeable reflections and shadows.

As stated before, medium and high powers have a comparatively short working-distance, and the front lenses will be so close to the water, urine, blood, etc., that the capillary attraction will often cause an adherence to the front surface of the objective; besides this, there is such a considerable depth to the fluid that it obstructs the light, requires a great change in adjustment for the various planes, and is usually in such vibration that a sharp focus becomes impossible; by merely dropping a cover-glass upon it all these objections are overcome.

If, however, this does not give sufficient space between the slide and cover, as, for instance, in the observation of living forms, a *live cage*—(Fig. 22)—or *compressor*—(Fig. 23)—should

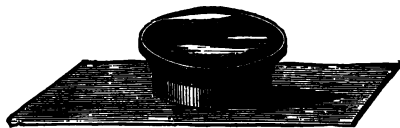


Fig. 22.

be used. In either of these the space may be increased at will; in the first by a sliding arrangement, in the second by means of a screw.



Fig. 23.

The above are merely practical considerations, but there are others of a theoretical nature and of as much importance. After a high power objective has been corrected to a certain thickness of cover, any variation, not necessarily considerable, has an injurious effect upon the spherical corrections, and consequently upon the resolving power. It is manifest that the quality of the latter will decrease as the variation increases, and when it reaches a point where no cover is used, it may be so considerable as to destroy an accurate perception of what is sought.

In this connection it is considered important to state what thickness of cover-glass it is best to use. As is probably well known, there are three grades, which are designated as No. 1, No. 2 and No. 3. Although they are classified, those of the same number are not absolutely of one thickness. The variation is about as follows: No. 1, 1-150 — 1-200 inch thick; No. 2, 1-100 — 1-150 inch thick; No. 3, 1 50 — 1-100 inch thick. According to the

prices of cover-glasses, when purchased by weight, the No. 1 give the greatest number and No. 3 the least. It may for this reason be thought, that the purchase of No. 1 is most advantageous, but it must be considered that there is a greater proportion of breakage by cleaning, as they are very thin and sensitive. Considered only from an optical standpoint, No 2 should generally be used, as the medium and high power objectives are adjusted to this thickness and give the best results with the thinnest of these. The same thickness is also used on test objects, but they are generally not of as much uniformity, as might be desired. Objectives sometimes have such an extremely short working-distance, that it is necessary to use the thinnest of No. 1, but as these are usually provided with adjustment for correction, their injurious influence is not so much felt. The thickest covers are most comfortable to handle and may be used with low power objectives without much sacrifice of definition.

To Draw Objects.—It is very important that the appearance of an object should be put upon paper, especially of one which is not permanently mounted. To do this does not require any great amount of skill, as the lines which are projected upon paper are merely followed out; but it is necessary that those drawings be made truthful. Nothing should be put down which is not actually seen;

neither should anything be omitted. Drawings thus made form a valuable record, not only for the individual, but for others who are following the same line of study.

Camera Lucida—How to Use It.—The camera lucida is the apparatus by means of which drawings are made. There are various forms of these, but those in most general use are the Wollaston Prism and neutral tint glass. The methods of attaching to the tubes are also numerous, but a very simple and effective device is that shown in Fig. 24. The



Fig. 24.

mounting is made of hard rubber, and is fixed to the cap of the eye-piece by means of a flexible, grooved ring.

The procedure of working should be about as follows: Focus upon the object and then incline the body, so that the center of the eye-lens will be 10 inches from the table. To obviate repeated measurements, a standard stick of this length may be used. If the instrument is so low that it will not allow the inclination of the body to an angle of at least 45 degrees when at this distance, it should be placed upon a box; or, if not too high, upon the case of the microscope. Now readjust the mirror and attach the camera lucida from below and place the paper under the instrument; look into the camera lucida from above, being careful that the eye is directly over the center of its opening, and the image will be found to be projected upon the paper. Possibly, and very probably, it will appear faint. This is due to the fact that the paper is almost as highly illuminated as the field. To remedy this defect a cardboard should be placed between the paper and light, so that the former will be shaded; the object will now come out in strong contrast. Take a well pointed pencil and follow the lines in the image. A little difficulty may at first be experienced in seeing both the pencil and image at one time, but after a little practice this is overcome.

It very often occurs that the pencil point can on no condition be seen distinctly, but this is usually due to abnormal sight, with which persons are often afflicted. In these cases, the glasses which are required for reading should also be used in drawing. The difficulty is not experienced in the image, as this can be adjusted to the eye.

Determining the Magnifying Power.—Although a magnifying table may be furnished, this gives the powers merely approximately, as more or less variation occurs in objectives and eye pieces of the same kind. As it is interesting, and sometimes important, to know the exact magnifying power, a simple method is mentioned. Procure a *stage micrometer*, divided into 1-100 inch and 1-1,000 inch, and perhaps 1-5,000 inch, or, if preferred, any suitable division in millimetre. Place the micrometer on the stage, focus and incline the microscope, as if for drawing, to within 10 inches from the table and attach the camera lucida; for low powers, 1-100 divisions may be used; for the higher ones, 1-1,000 or more; the divisions as now projected may be marked upon the paper and then measured off with a rule divided into inches and 1-10 inches; if, for instance, the 1-1,000 divisions are used, and one division on paper covered 1 inch on the rule, it is evident that the magnifying power is 1,000 times; if it covered 2-10 (=1.5) on the rule, it would be 200 diameters, and so on.

Measuring the Size of an Object.—A simple and reliable way of learning the size of an object is by means of an eye-piece micrometer. As, however, this does not measure the object directly, but only its image, the first part of the process makes it more complicated. However, this portion is usually attended to by the manufacturers. The eye-piece is provided with a slot, into which a micrometer is fitted. A micrometer with the same divisions as the eye-piece micrometer is placed upon the stage and the objective focussed upon it. It is now observed how many of the divisions of the eye-piece micrometer are contained in the magnified division of the stage micrometer, and the resulting figure is placed in the sub-division under the objective. To determine the actual size of an object, this is now placed on the stage and, noting the number of divisions which cover it, these are divided by the number on the card, and the resulting figure gives the actual size. Suppose the figure on the card is 8.0 and the image of the object covers 40 of the spaces which are divided into 1-10 milli-metre, the size of the object would be 5-10, or 1-2 milli-metre; or, expressed in inches, (25.4 m. m. = 1 inch) about 1-50.

ADVANCED MANIPULATION.

Dry Adjustable Objectives.—The information which has thus far been given on the manipulation of the microscope may be termed initiatory, as it is supposed (at any rate hoped) to have disclosed some new principles. These are comparatively simple, and, with a moderate amount of attention, are easily acquired. It is the intention now to speak of something more complex and to give instruction in the use of higher grade adjustable and immersion objectives. The difficulty of doing this increases with each step of advance, and whether it can be overcome by means of written words, is, perhaps, an open question. However, the writer is certain that if the following instructions are faithfully adhered to, satisfactory results will be gained. The highest attainment must of necessity be the result of perseverance and knowledge of the various properties of an objective, which are given in preceding pages.

It is assumed that a dry objective is used, say a 1-6 140° or 1-8 135°, and provided with screw-collar adjustment, for cover-glass thickness; it is further assumed that in the first, the variation between the two first systems (anterior and middle) is attained

by means of a rectilinear motion to the middle and posterior system and stationary anterior system, while in the second the conditions are reversed; the two posterior systems remain stationary while the front is adjusted. Both are arranged with graduations upon the screw-collar and have an index. In both zero represents the adjustments as *open, i. e.*, the lenses are separated to their fullest extent, and the objective is approximately adjusted for the *thinnest* covers. As the adjustment is moved towards the higher numbers, it is *closed*, and indicates correction for gradually thicker covers. At the medium point 5, the correction is about for medium thickness, while at 9 it is for the thickest. Before the objective is attached the adjustment should be closed, as, if this is neglected and the objective has a short working distance, the front lens may come in contact with the cover when it is focused.

Probably the best object for studying the effect of the screw-collar adjustment and acquiring skill in determining its best point is again *P. angulatum*. Place a slide of this upon the stage, and with a low power eye-piece select a diatom which appears to be flat; such a one may usually be found when there are a number on the slide. After the objective has been attached to the nose-piece, focus carefully and observe whether any lines can be seen; if not, grasp the *milled edge* of the *adjustment collar* between the thumb and first finger of the *left* hand, keeping the fingers of the *right* hand upon the micrometer screw,

or *vice versa*, if from the outset it was made a habit to use the left hand on the fine adjustment; turn the collar slightly toward its open point, and as this will place the object out of focus, move the fine adjustment correspondingly; continue to turn the collar, little by little, and do not cease to observe closely; also, after each movement, focus above or below the plane of the object, so that this will be indistinct, and look for the lines. Possibly after a little they will begin to appear faintly; but, if not, continue to bring the collar toward the middle. The lines must now soon make their appearance, and, when they do, it will probably be above the plane of the diatom. This is an indication that the objective is approaching its correction for the cover. Now keep the *lines* in focus, while the correction collar is being gradually turned, until the lines and the outline of the diatom lie in one plane; the objective is now said to be *corrected for cover*. Observe which number corresponds to the index, and note this upon paper; again return the collar to its closed point and go through the same proceeding as carefully as at first. When the best point is again reached, look for the number and see whether it agrees with the first; very likely it does not, which is owing to a want in the faculty of perception, due to a too slight acquaintance with the phenomena. These trials should be repeated until the proper sensitiveness of feeling in making the adjustments is acquired, and until they can be made to correspond with a certainty to within at least two divisions.

Remove the eye-piece and attach one of higher power. It must now, however, be remembered that if there is a considerable difference in the powers, there will be a relative difference in their lengths, and that this will cause a difference in the optical length of the tube; this not only will require another adjustment for focus, but will partially destroy the correction as made with the low power. After some practice, the amount of variation may be fixed upon and may be noted for the future; but, to determine it, the same plan as suggested with the low power eye-piece should be followed.

When it is found after repeated trials that sufficient skill has been acquired to bring the collar to within one division, the number and power of the eye-piece should be scratched with a diamond upon the slide or with pen and ink upon the label; thus, if it is found that with a $1\frac{1}{2}$ inch eye-piece the index shows 5, and with a $\frac{1}{2}$ inch eye-piece shows $5\frac{1}{2}$, it should be marked $1\frac{1}{2}$ —5 and $\frac{1}{2}$ — $5\frac{1}{2}$. For future examinations on the same slide, this will facilitate work and give the assurance that the best results are thus gained without further trial.

Mr. Wenham's general rule for obtaining the best correction on objects in general is as follows: "Select any dark speck or opaque portion of the object and bring the outline in perfect focus; then lay the finger on the milled head of the fine adjustment and move it briskly back-

wards and forwards in both directions from the first position. Observe the expansion of the dark outline of the object, both when within and when without the focus. If the greater expansion or coma is when the object is without the focus or farthest from the objective, the lenses must be placed farther assunder (or opened). If the greater coma is when the object is within the focus, or nearest the objective, the lenses must be brought closer together (or closed). When the objective is in proper adjustment, the expansion of the outline is the same both within and without the focus."

Immersion-Adjustable Objectives.—As was stated before, immersion contact between the objective and cover glass is made by either water or homogeneous fluid. The fluid should be kept in a small bottle or phial, the cork of which is pierced to receive a small pointed stick or match, and this should project sufficiently so that it will enter the fluid about $\frac{1}{2}$ -inch. The fluid will then always be free from dust, and by withdrawing the cork the stick will always carry a drop of fluid with it.

In fixing an immersion objective to the stand, the latter should first be put in an upright position; the fluid should now be attached to the front lens, but care should be taken not to put on too much; it should be merely enough to cover the surface. If too much, portion of it should be removed by

allowing it to adhere to the finger. The objective may then be attached to the stand and brought down until the fluid is in contact with the cover; the stand is now inclined and the objective focused; if this method is followed there is no danger of flooding the entire cover with fluid, which sometimes may be the means of destroying the object; neither can the fluid run out from between the two surfaces.

Extreme cleanliness should be observed in all work connected with the microscope, and particularly in the use of immersion objectives. The use of immersion fluid in itself involves a certain amount of inconvenience, but as in many cases it is absolutely necessary, the observance of fixed rules will materially help to overcome some of the disagreeable features. After the work with an immersion objective has been completed, the objective should be removed from the stand and its front, as well as the slide, should *invariably be cleaned*; the fluid may be removed by a moist piece of soft linen and then cleaned with a dry piece; chamois skin is not suitable, as it does not absorb the fluid.

Test-Plate.—Almost all microscopists who take an active interest in the capacity of their instruments, supply themselves with a set of *test objects*, of which *P. angulatum* is in most general use or with a so-called test-plate. These plates consist either of a series of bands of finely ruled lines ranging from

5,000 to the inch to 120,000 to the inch (beautiful specimens of these are made by Prof. Wm. H. Rogers, of Cambridge, Mass., and C. Fasoldt of Albany, N. Y.,) or with a series of diatoms, upon which the markings represent certain divisions of an inch. The one of these which is principally used is made by J. D. Moeller and consists of a series of 20 diatoms. They are furnished mounted both *dry* and *in balsam*, but the latter is the most common. Below is a table giving the numbers, names of the various diatoms and divisions on their surfaces to 1-1,000 inch. A specimen of *Eupodiscus Argus* begins and ends the series:

Striæ in 1-1000
of an inch.

| | |
|--------------------------------------|--------------|
| 1. Triceratium Favus Ehrbg..... | 3.1 to 4. |
| 2. Pinnularia nobilis Ehrbg | 11.7 to 14. |
| 3. Navicula Lyra Ehrbg. var | 14.5 to 18. |
| 4. Navicula Lyra Ehrbg..... | 23. to 30.5 |
| 5. Pinnularia interrupta Sm. var.... | 25.5 to 29.5 |
| 6. Stauroneis Phoenicenteron Ehrbg | 31. to 36.5 |
| 7. Grammatophora marina Sm | 36. to 39. |
| 8. Pleurosigma balticum Sm | 32. to 37. |
| 9. Pleurosigma acuminatum (Kg.) | |
| Grun | 41. to 46.5 |
| 10. Nitzschia Amphioxys Sm..... | 43. to 49. |
| 11. Pleurosigma angulatum Sm..... | 44. to 49. |

| | | |
|----------------------------------------------------------------------------------------|-----|--------|
| 12. <i>Grammatophora oceanica</i> Ehrbg. =G. <i>subtilissima</i> | 60. | to 67. |
| 13. <i>Surirella Gemma</i> Ehrbg..... | 43. | to 54. |
| 14. <i>Nitzschia sigmoidea</i> Sm..... | 61. | to 64. |
| 15. <i>Pleurosigma Fasciola</i> Sm. var... | 55. | to 58. |
| 16. <i>Surirella Gemma</i> Ehrbg..... | 64. | to 69. |
| 17. <i>Cymatopleura elliptica</i> Breb..... | 55. | to 81. |
| 18. <i>Navicula crassinervis</i> Breb.= <i>Frus-</i> <i>tulia saxonica</i> Rabh.... | 78. | to 87. |
| 19. <i>Nitzschia curvula</i> Sm. | 83. | to 90. |
| 20. <i>Amphipleura pellucida</i> Kg..... | 92. | to 95. |

It may be said, and, perhaps with truth, that a test plate does not belong to the necessities of an outfit, but considering that it is a guage, on which the optician usually bases the quality of his objectives, it is valuable to the owner of an objective to be able to determine whether, under his manipulation, the objective will perform as well as is claimed for it; due consideration must, however, be given to the fact that there is a certain amount of variation among different plates, as is shown in the above table. Outside of this, it is a continual incentive to determine the extreme performance of an objective and it thus becomes the means of acquiring great manipulative skill, which can not be underrated. The writer is in a position to know that there is great need of this; innumerable cases have come to his notice where several objectives

of the same kind and equal quality gave unequal results in different hands, and would be highly eulogized by the possessor of one and condemned by that of another.

Immersion Objectives on Test Plate.—To determine the highest capacity on test objects, ordinary day light is hardly sufficient; moderate sun light or good lamp light is best suited, but the latter, from the fact that it is always at hand, is preferable. For the purpose of explanation, we will assume that a flat-wick lamp and a 1-8 or 1-10 homogeneous immersion objective is used. If the right hand is used on the micrometer screw, place the lamp at the right side of the instrument, about 10 inches from it, with the *edge* of the flame turned toward the mirror.

The test plate may now be placed upon the stage, and as the diatoms in balsam are very transparent, and therefore very difficult to find, a lower power objective may be used as a finder; bring No. 1, or *Triceratium Favus*, in the center of the field, and after the objective has been removed, attach the immersion objective in the manner prescribed; the adjustment collar may be placed at zero, as this is about the correct point for standard length of tube. Get the best possible illumination with the mirror at the central point and move the test plate from diatom to diatom until it reaches No. 11, *P. angulatum*, but observe closely the structure of each one as it comes into the field. Next see whether the objective is corrected; if the lines and outline, or

middle rib, do not appear to be in one plane, adjust the collar until they are, and then continue the advance toward the higher numbers until one is reached on which no lines can be seen. Swing the mirror-bar to an obliquity of 20° . and, readjusting the mirror, observe the effect. It is very probable that the lines will show, and, if so, continue the advance; if they do not, give 10° or 20° more obliquity, and after the structure comes out, again go forward. A point may thus be reached, where with the greatest obliquity which can be given and with the best possible illumination, the objective seems to have come to the limit of its performance. From the claims which have been made for it, it ought to do better. What is the cause of failure? Possibly the mirror is not correctly focused, or the adjustment collar may not be correct for oblique light; perhaps the eye-piece does not give sufficient magnifying power to distinguish the striæ. It may be any one of these causes or all combined. As to the eye-piece, the manipulator must remember the amount of separation of lines in the last object which was resolved, and from the gradation in the coarser specimens must judge whether the power is sufficient; it should be added that for any over No. 14 and under No. 18 a $\frac{3}{4}$ inch eye-piece should be used, and for those above No. 18 a power of $\frac{1}{2}$ inch will probably be necessary, provided a $\frac{1}{8}$ or one-tenth objective is used. After this condition has been complied with, look to the correction col-

lar of the objective ; to obtain the highest results it very often occurs that a different adjustment is required for oblique light from that for central light. Note the number at which it stands, and then work it back and forth, watching carefully for results. If this has no influence, return it to its number or to a point where the outline of the object appears most sharp. Now look to the illumination ; vary the distance of the mirror to the object, and, if it conflicts with the stage or does not give the desired results, vary the distance of the lamp to the instrument and watch the effect of the change through the tube. A great change in the illuminating power can thus be produced ; the light is best when it covers the least space, as it is then most intense. The light may be quickly adjusted by throwing it upon a point on the slide in the opening of the stage and watching it there. If neither of these changes give any improvement, recourse must be had to another expedient. Place a bulls-eye between the lamp and the mirror in such a position that the light will be well thrown upon the latter. For this purpose it should be moved back and forth. Keep it a little below the line of the face of the stage, so that the light will not strike it on its upper and as little as possible on its lower surface ; if the light from the bulls-eye directly reaches the object, it destroys the effect of the oblique illumination. Great care should be given to this point, as it is very important.

If all of these suggestions have been followed, a great difference will undoubtedly be noticed in the performance of the objective ; but if it still does not come up to the standard, patience must not be lost. The slightest change in the mirror, bulls-eye, or lamp, a touch to the correction collar or micrometer screw is sometimes followed by astonishing results. The beginner should sit down with the expectation that he will fail at the first trial. At each succeeding trial he can easily notice his improvement in manipulation and the gain of corresponding results. He should be able to bring the performance of the objective up to the claims made for it, if it has come from the hands of a reliable optician, and should not rest until this is accomplished.

To the histologist it may seem strange that the writer has thus far only spoken of working with objectives on diatoms. This, however, was done advisedly. They are thin, and therefore as suitable as a thin section and far more preferable than a thick one. Their form and structure is easily recognisable, and there is very little variation among those of the same kind ; therefore, rules laid down regarding them are generally good. It is conceded by advanced workers that the time spent over diatoms for the purpose of studying objectives is well applied, and the most expert manipulators have acquired their experience in this manner. An objec-

tive which works well on diatoms works equally well on other objects, and therefore the manipulative skill which has been attained on the former is as well applied on the latter. At the outset, work may be done on other objects than diatoms, and where ordinary working objectives, such as a student 1 inch and 1-4 inch, or 3-4 and 1-5 inch comprise the outfit, the road to good manipulation may be as short as with diatoms. The conditions in both cases remain the same; but it must be cautioned that, if histological preparations be used, only such be selected as are reliable. A poor specimen is perhaps as bad as none at all; an abnormally thick one obstructs light, makes it impossible for the objective to penetrate through the various layers, and leaves the impression that the latter is defective.

SUB-STAGE ILLUMINATION.

This chapter should form an important part in a book of instruction on the use of microscope. For many years after the discovery of the achromatic objective, sub stage illuminators formed an important part of an instrument, as the loss of light with medium and high powers was considerably greater than it is to-day; also because the method of mounting the mirror was more primitive than it is at present. Then after new methods were devised for increasing the angular aperture and thus attaining a higher performance in objectives, the necessity for it was not so much felt, and it came into disuse to a certain extent. Of late years, however, the demand for it has again sprung up, but in a form to meet present requirements.

By means of a condenser of proper construction, better results are reached than by the mirror only, and although a certain amount of delicacy of manipulation is necessary, the best performance of the objective is accomplished more easily with one than without one. Very often results may be gained with one of good construction, which it is impossible to reach without one.

Objectives and Eye Pieces.—A medium power objective does good service as a sub-stage condenser, but its length precludes its general use. To attach it to the microscope it should be screwed into a sub-stage adapter with society screw and fitted to the sub-stage with its front system toward the object.

The Periscopic eye-piece is similar in its construction to the so-called Webster condenser, and gives excellent results. It is attached in the same manner as an objective, except that a special adapter must be provided. In the Professional microscope of the Bausch & Lomb Opt. Co., the diameter of tube corresponds to that of the sub-stage, so that the eye-pieces can be attached without special fitting.

Hemispherical Lens.—This is probably the most simple form of sub-stage illuminator. It is a plano-convex lens of about $\frac{1}{2}$ -inch diameter, and is attached directly to the slide by means of an immersion fluid. Where no great obliquity is employed, a better plan is to cement it to a thin slide and then place the slide with object upon this. It is not as its shape might imply, a condenser, except when the plane mirror is used, when it partially acts as such; theoretically the light from the concave mirror should pass through it without undergoing any change. A greater obliquity of light is obtained by it, than is possible without it.

Wenham Button.—This is the invention of Mr. Wenham, and is attached to the slide in the same manner as the hemispherical lens. It is a semi-circular glass with two flat surfaces; the third and upper surfaces are polished and the circular one is ground and polished on a curve of a shorter radius than the circle, thus permitting a certain amount of concentration of light.

Woodward Prism.—This device was suggested by the late Col. J. J. Woodward, and is attached as mentioned above. It is a triangular prism and is intended to allow the use of extreme oblique light. It is in no sense a condenser, but merely permits oblique light to reach the object without undergoing refraction.

Narrow Angle Condensers.—These may be divided into two classes—achromatic and non-achromatic. In effect they are similar to objectives, but are better adapted to the purpose, from the fact that they are provided either with revolving diaphragms or with caps to decrease the aperture and thus the amount of light. The most simple forms are of such a length that they can be used on instruments, in which the sub-stage is fixed to the stage, but the more complex and all those of which descriptions follow, should be used on adjustable sub-stages.

Wide Angle Condensers.—The most simple form of these is the so-called *Immersion illuminator*.



Although it is intended to be used in connection with the slide by fluid it may also be used dry, as may also all immersion condensers; but it must be remembered that it as well as others give a greater amount of light when used with immersion fluid, as all rays reach the object, which is not the case when used



Fig. 25.

dry. It may be used in plain adapter, when all the light which reaches the lower lens (which is nearly the diameter of the Society screw) is concentrated to one point. Another mounting is provided, by means of which any amount of oblique light may be obtained, without changing the position of the mirror. This feature is particularly valuable in connection with those stands in which an extremely thick stage or the construction of the mirror preclude any considerable obliquity of light.

Ward's General Illuminator.—This primarily is similar to the above. The same optical portion is used, but it is mounted in a sub-stage fitting. To the lower side of this an Iris diaphragm is fitted and made to traverse by means of a screw, from the extreme edge to the center. The amount of illumination

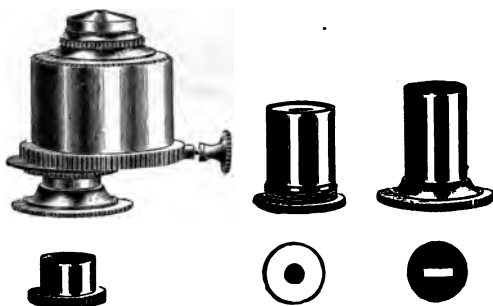


Fig. 26.

and obliquity is thus under constant control. The Iris diaphragm may be removed and in its place may be attached other adapters with various stops.

Large Sub-Stage Condenser.—The construction of the optical portion of this accessory is similar to that of the Abbe condenser, with which Dr. Koch made his original investigations, and in its results is identical with it. It is the most approved form at the present time, not only on account of its wide angular aperture, but on account of the extremely large size of the lenses, which receive all the light which may pass through the sub-stage. It may be used dry, but is intended to be in homogeneous immersion contact with the lower surfaces of the slide. Three styles of mounting are provided, one with plain adapter, another with a mounting—(No. 1, Fig. 27),—which, by revolving a milled ring, causes a diametrical movement of a $\frac{1}{4}$ inch aperture from the extreme edge



No. 1.



Fig. 27.

No. 2.

to the center, thus giving all grades of obliquity from zero to the highest point. The third—(No. 2, Fig. 27)—is provided with a swinging ring, which in its recess receives diaphragms of various size apertures. Some of these decrease the volume of light; others with central stops give dark ground illumination; another has an eccentric aperture, by means of which various grades of obliquity of light may be attained. Besides these a blue glass is provided for monochromatic light with lamp light, which may be used alone or in combination with any of the diaphragms.

This condenser is one with which the highest results are gained, and which combines in itself the means of accomplishing almost all the various illuminating effects.

In using any of the condensers which are attached to the sub-stage, different results are gained by varying their distance to the object. This becomes necessary on account of the difference in thickness in slides, to properly focus the condenser; but in many cases the best results are reached when the light is not intense. A proper judgment of the best point can only be reached after repeated trials.



CARE OF A MICROSCOPE.

Besides acquiring the ability to properly use an instrument with its accessories, it is important to know how to keep it in the best working condition. It may be said without reserve that an instrument properly made at the outset and judiciously used should hardly show any signs of wear either in appearance or in its working parts, even after the most protracted use; and further than this, every good instrument should have a provision for taking up lost motion, if there is a likelihood that this may occur in any of the parts.

Especial care should be given to the optical parts, in fact such care, that they will remain in as good condition as when first received. Accidental injury may of course occur to them, but if a systematic manner of working is followed and a special receptacle for each part is provided, this may usually be avoided. The following rules refer mainly to the instruments manufactured by the Bausch & Lomb Opt. Co., but are applicable to instruments in general.

To Take Care of a Stand.—One of the first rules should be to keep the instrument *free from*

dust. This may be done in a manner formerly prescribed. If dust settles on any part of the instrument, remove it first with a camel's hair brush, and then wipe carefully with a chamois skin, *with* the grain of the finish and not across it, as in the latter case it is likely to cause scratches. Keep the working and sliding parts absolutely free from dust, as this grinds and will thus soon cause play.

Use no alcohol on any part of the instrument, as it will remove the lacquer. As the latter is for the purpose of preventing oxydation of the metals, it is important to observe this rule.

In using draw-tube, impart a straight up and down motion to it. If a spiral motion is given to it, it will cause scratches, and in time will wear off the nickel.

If it becomes necessary to lubricate any of the parts, use a slight quantity of soft tallow or good clock oil.

In an instrument which is in constant use, it sometimes occurs that the pinion works loose and occasionally to such an extent that the body drops of its own weight. Tightening screws are provided to take up the play—in the Professional, American Concentric, Universal, Griffith Club and Physician Microscopes these are at the back of the pinion. In the Investigator, Model and Family Microscopes, they are seen in the slide by removing the body.

In using a screw-driver, grind its two large surfaces so that they are parallel and not wedge-shape, and so it will exactly fit in the slot of the screw-head.

In inclining the stand *always* grasp it at the arm and *never* at the tube, as in the latter case it may loosen the slide or tear off some of the parts.

When repairs or alterations are necessary, always have these made by the manufacturers; they can, from the system of duplicated parts, not only do it cheapest, but best.

To Take Care of Objectives and Eye-Pieces.—

It is as necessary to keep these *free from dust* as the stand, in fact even greater cleanliness should be observed. When indistinct, dark specks show in the field, the cause may usually be looked for in the field-lens, although sometimes in the eye-lens also. The dust may be removed by a camel's-hair brush, but when this is not sufficient, use a well washed piece of linen, such as an old handkerchief. From its fine texture chamois skin is desirable, but as it is fatty it should never be used until after it has been well washed.

The same method applies to cleaning objectives. Clean an immersion objective invariably after it has been used, first by removing the fluid by a moist linen and then by using a dry piece.

Keep the objectives especially in such a place where they are not subject to extreme and sudden changes of temperature, as the unequal expansion and contraction of glass and metal may cause the cement between the lenses to crack. Also keep them from direct sun-light.

Screw them into the nose-piece and unscrew, by grasping the milled edges.

Avoid any violent contact of the front lens with the cover-glass. Usually the latter suffers, but it is as liable to occur to the former.

Above all, it should be made a rule that no one but the owner handle the microscope and accessories. One person may be expert in the manipulation of one instrument and still find it difficult to work with another. The fine adjustment particularly causes the greatest difficulty as in some instruments, it corresponds with the movement of the micrometer screw, while in others it is contrary and thus the objective as well as object are endangered.

APPENDIX.

(Re-printed from "The Microscope," Jan., 1885.)

CONSIDERATIONS IN TESTING OBJECTIVES.

EDW. BAUSCH.

There is a laudable desire in almost all persons possessing a microscope to become intimately acquainted with it, and for this purpose it is not only necessary to learn the use of its mechanical parts, but to understand its optical capacity, which is considerably more difficult, and which involves more considerations than would appear on first thought.

With all the care which may be bestowed upon objectives, they are to a certain extent, works of chance, and depend upon the optician's judgment, industry and skill, and upon the variations in glass, for their excellence and uniformity. These conditions are often so varying that in the case of sev-

eral objectives of the same formula, made at the same time there will be such great differences that it can hardly be conceived on the first examination, that they were to be similar. It is at this point especially necessary to detect the errors, to determine their cause and apply the remedy, and do this properly often involves an inconceivable amount of work, and in many cases the final results are reached at a pecuniary loss. There are certain fixed tests for each kind of objective, and to the best of my knowledge all reputable opticians bring each objective up to its standard before allowing it to pass their hands, irrespective of the cost of doing so. This must of necessity be so, if only out of business consideration, and not for a love of each production, for it is evident that a well-earned reputation would soon lose its pre-eminence, and would acquire one for unreliable or poor work, if on comparison, objectives of the same kind would show a marked difference. There is sometimes a fortunate combination of circumstances which makes a certain objective better than its fellows, but this is a rare exception, and is positive evidence that the exact requirements of the formula have been complied with. As a rule, therefore, I believe that the opticians' claim may be relied upon, and where the results in the hands of the microscopist do not correspond with them, the cause may usually be looked for in the lack of experience in manipula-

tion or in conditions, which differ from those under which the objective was completed. The belief, which I am aware is extant, that there are great differences in objectives purporting to be similar, is, in my opinion, not justified, at any rate in the productions of those men who, by general acknowledgement, are at the head of their profession. I admit that, as in everything which depends upon human skill, there is strictly speaking, no absolute uniformity, but also claim, that with few exceptions, the differences are so slight, that anything but the most expert manipulation cannot detect them.

It therefore appears to the writer that any information which will tend to improve the knowledge of testing objectives will not only prove beneficial to the microscopist, but will prove advantageous to the optician, in that his work will receive a fair trial, based upon a knowledge of the principles involved, and that he may be convinced that all his work which deserves commendation will be the better appreciated. The following points are by no means new, but are often lost sight of in making tests. The writer will speak of medium and high power objectives only, as the deleterious influences are most noticeable in these, but they apply as well to the lower powers though in a less degree.

The part of the instrument which has a strong bearing in the performance of the objective is the mirror. It should be adjustable on the mirror bar,

so that it can be accommodated to the variations in distance of the source of light from the instrument. When parallel rays are used, as with light from the sun or clouds, its distance from the object should be decreased and increased when lamp-light is used. It should be exacted that the focus of the concave mirror be within the limits of its adjustment. The serious disadvantage of its incorrectness in this respect can easily be seen by taking, for instance, a 1-5 objective which will resolve *P. angulatum* nicely with central light, when the mirror is exactly focused. By moving the latter out of focus it will be seen that the objective loses in performance, and if this is carried sufficiently far it will arrive at a point where the objective will cease to show any lines. The effect will be the same on any other object, and is caused by the lack of proper concentration of light on the slide. When oblique light is used, unless the diaphragm moves with the mirror, it should be removed, as the advantage of obliquity is diminished or destroyed by the loss of light.

The cover glass exerts probably the greatest influence in testing as well as in general work. This should be used of a thickness which corresponds to that to which the objective (if non adjustable) was originally corrected. If thicker or thinner covers be used, the objective will be spherically over or under corrected, and will have to be moved corresponding above or below the plane (outline) of the ob-

ject to distinguish its structure, if the variation is considerable the difference between the two planes will be so great that it will cease to show any structure, and it may then be said to be lacking in defining power although in reality it possesses it, but is not properly used. Generally speaking the objective may be said to be spherically corrected when it gives the best defined image; that is, when the outline and internal structure of an object of extreme thinness appear in one plane. When, after the objective has been focused on the outline of the object, it is necessary to increase the distance to focus on the structure, it is evidence that the objective is spherically over-corrected and that the cover is too thick; in adjustable objectives the correction collar must be brought to its closing point, which means that the lenses are brought in closer contact. When the objective must be focused to a point beyond the outline of the object to see its structure; that is, brought closer to the cover-glass, it proves that this is too thin, and is then said to be spherically under-corrected; to give the proper adjustment in an adjustable objective in this case the adjustment is opened—the lenses are separated. It requires a certain amount of study to distinguish these phenomena, and although it can be done in well prepared specimens, I know of none better than coarsely marked diatoms, such as *P. angulatum*.

Although I am aware that many eminent microscopists do not favor adjustable objectives for every

day work, I must confess that I fail to see the force of their arguments. From the foregoing it will be seen that unless the cover-glasses are of a thickness corresponding to that which was originally used, the objective may be made to do imperfectly what is in its power to do well, and when pressed to its full capacity may and is likely to fail. It must be remembered that cover-glasses of the same number are not of the same thickness. The selection of those of proper thickness is expensive and tedious, whereas the knowledge of correcting the objective is easily acquired, and in the latter case it is in the manipulator's power to command the highest performance of which the objective is capable; further than this, it has the advantage that it may be used as a non-adjustable objective if desired. When homogeneous immersion objectives were first introduced they were mounted in fixed settings, as it was expected that the thickness of the cover-glass would not affect the correction; although this assumption was correct, it was found that even in these it was necessary. How much more then, is it required in dry or water immersion objectives?

Another factor in the disturbing influences is the variation in length of tube; the deleterious results are similar to those with varying cover-glasses. Objectives are usually adjusted to $8\frac{1}{2}$ or 9 inches length of tube, and although this in itself is a fixed standard, it usually becomes variable by changing objectives and eye-pieces. That this is so in objec-

tives is patent, and that it is so in eye-pieces can easily be determined by making a change in powers, when it will be found that a change in focus is required. By decreasing the length of tube the objective will appear to be spherically under-corrected and *vice versa* when it is increased, so that it is apparent that by the use of the draw-tube the effect of the cover-glass may be partially neutralized; for instance, when by the use of a thin cover the objective is spherically under-corrected, it may, to a certain extent, be corrected by causing a corresponding over-correction in the tube by increasing its length. The use of the draw-tube for the purpose of changing the amplification or for the matter of convenience can hardly be commended, except in cases where adjustable objectives are used.

Considerable also depends upon the perfection of the eye-piece. I believe that, as a rule, too little care is devoted to it; at any rate, it is certain that while any Huyghenian eye-piece for a telescope can be used on a microscope, very few which have been made for this can be used on a telescope; and while it is true that no such perfection may be required in the former, it leaves such an indefinite range that it may become difficult to place a limit for the perfect and imperfect. In all work, and especially in testing, it should be seen that the eye-lens, as well as the field-lens, are perfectly clean,

Among the absolutely necessary conditions in judging of the quality of an objective are perfect specimens, especially if they are sections. A thick object obstructs the light and generally makes it necessary to go through so many layers or planes that it is difficult to get any one distinct ; the impression may thus easily be given that the objective is at fault. The difference between two objects of the same nature may be so great that, while with one the objective may be condemned as imperfect, it may with the other appear to be of extraordinary excellence.

In conclusion, I will say that there may be other conditions which may influence the performance of a lens, and to acquire the power of eliminating them requires considerable experience. When an objective does not correspond with the claims of the optician, judgment should not be passed upon it until after repeated trials have been made, in all of which the above points should not be lost sight of.

THE
BAUSCH & LOMB OPTICAL CO.

HAVE THE MOST

IMPROVED MECHANICAL FACILITIES,

AND ARE THEREFORE ABLE TO MAKE

ANY KIND OF

Special Apparatus

AND

Mechanical Devices

IN ANY DESIRED QUANTITY.

They will furnish ESTIMATES on receipt of
Drawings and Descriptions.

*On the following pages will be found
Illustrations of a number of Instru-
ments, manufactured by the BAUSCH
& LOMB OPTICAL CO.*

*For complete descriptions of Micro-
scopes; as well as OBJECTIVES and
ACCESSORIES, send for complete cata-
logue.*

Address :

P. O. Box 354, Rochester, N. Y.

P. O. Box 432, New York.



EXCELSIOR DISSECTING MICROSCOPE.



COMPACT DISSECTING MICROSCOPE.



FAMILY MICROSCOPE.




HARVARD MICROSCOPE.



MODEL MICROSCOPE.



INVESTIGATOR MICROSCOPE.

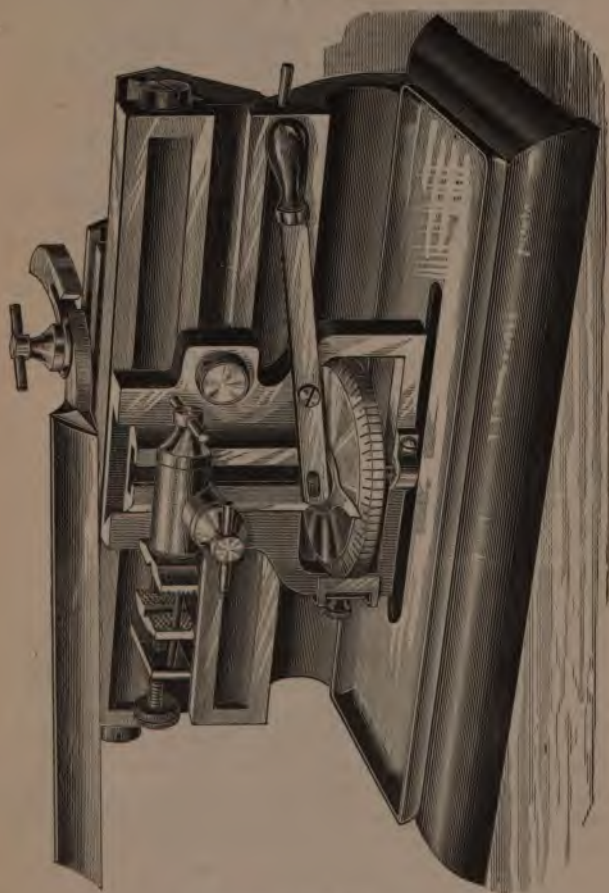

IMPROVED MICROTOMES.

THE MICROSCOPE

AN ILLUSTRATED MONTHLY JOURNAL.

EDITORS AND PUBLISHERS:

C. H. STOWELL, M. D., F. R. M. S.,

Prof. of Histology and Microscopy in the University of Michigan.

LOUISA REED STOWELL, M. S., F. R. M. S.,

Assistant in Microscopical Botany, University of Michigan.

Price \$1.00 a Year for the United States; 5 shillings for
Great Britain.

PUBLISHED AT ANN ARBOR, MICH.

"That Dr. and Mrs. Stowell have in them the elements necessary for successful journalism is apparent."—*MICH. MED. NEWS.*

"The names of the editors alone is ample guarantee to insure continued success to this journal."—*GOOD HEALTH.*

"It fills a sphere of its own and should be in the hands of every physician and druggist in the country."—*John Phil, in AM. JOUR. OF MIC.*

"It always has a neat appearance and is constantly filled with good matter."—*SANITARY NEWS.*

"It improves steadily and decidedly with every issue."—*BISTOURY.*

"It is produced in an excellent form, is well printed, and illustrated with good illustrations."—*SCIENCE.*

Twenty-two of the twenty-four pages of the January number for 1885, consisted of **ORIGINAL COMMUNICATIONS**; while the February number consisted **ENTIRELY OF ORIGINAL COMMUNICATIONS**, with the exception of two small items. In these two numbers alone there were **OVER TWENTY** original communications. The other numbers of the year have not appeared at this writing, but the peculiar feature of having the columns filled with the latest and best original matter will be continued through the year.

Address,

C. H. STOWELL,

ANN ARBOR, MICH.

AMERICAN MONTHLY MICROSCOPICAL JOURNAL.

ESTABLISHED 1880.

AN ILLUSTRATED PERIODICAL
FOR ALL WHO USE
THE MICROSCOPE

EITHER PROFESSIONALLY OR AS AMATEURS.

TESTIMONIALS.

Extracts from letters of well-known workers with the microscope. These being unsolicited and not intended for publication, the signatures are withheld.

"I hope the Journal will continue to prosper, and that its columns in the future, as in the past, will be filled with hints and suggestions to working microscopists, both amateur and professional."

"I enjoy your Journal very much. I learn many good things from it, and I am sure that is the case with all who read it."

"I want to thank you for your excellent Journal. I am pleased and satisfied and I want you to know it."

"I send you \$1.00 enclosed, and wish you much success for the new year and you richly deserve it."

"I am very well pleased with the Journal and its standing, considering the wide range of subjects it has to cover. But there is one feature about your Journal that cannot be praised too highly. * * * Although they might yield a nice income to the editor, I think this disinterestedness should be generally appreciated."

"I am glad to see by the last number that its present condition and future prospects are satisfactory; and for myself can say that it is always received and perused with pleasure."

SUBSCRIPTION PRICE, \$1 PER YEAR,
STRICTLY IN ADVANCE.

Editor and Publisher,

R. HITCHCOCK,

P. O. Box, 630.

Washington, D. C.

THE STUDENT'S MANUAL OF HISTOLOGY

FOR THE USE OF

STUDENTS, PRACTITIONERS AND MICROSCOPISTS

BY

CHAS. H. STOWELL, M. D.,

ANN ARBOR, MICH.

Price \$3.50.

PERSONS VISITING
NEW YORK CITY,

ARE CORDIALLY INVITED TO CALL AT OUR

BRANCH OFFICE,

37 Maiden Lane,

WHERE WE KEEP ON HAND


A COMPLETE STOCK

OF ALL ARTICLES MANUFACTURED BY US.


BAUSCH & LOMB OPTICAL CO.







[The page contains extremely faint and illegible text, likely bleed-through from the reverse side of the document. The text is scattered across the page and cannot be transcribed accurately.]



LANE MEDICAL LIBRARY

To avoid fine, this book should be returned
on or before the date last stamped below.

| | | |
|--------------------|--|--|
| DEC 24 1999 | | |
|--------------------|--|--|



C.M.C.